


Polarimeter tubes and Accessories

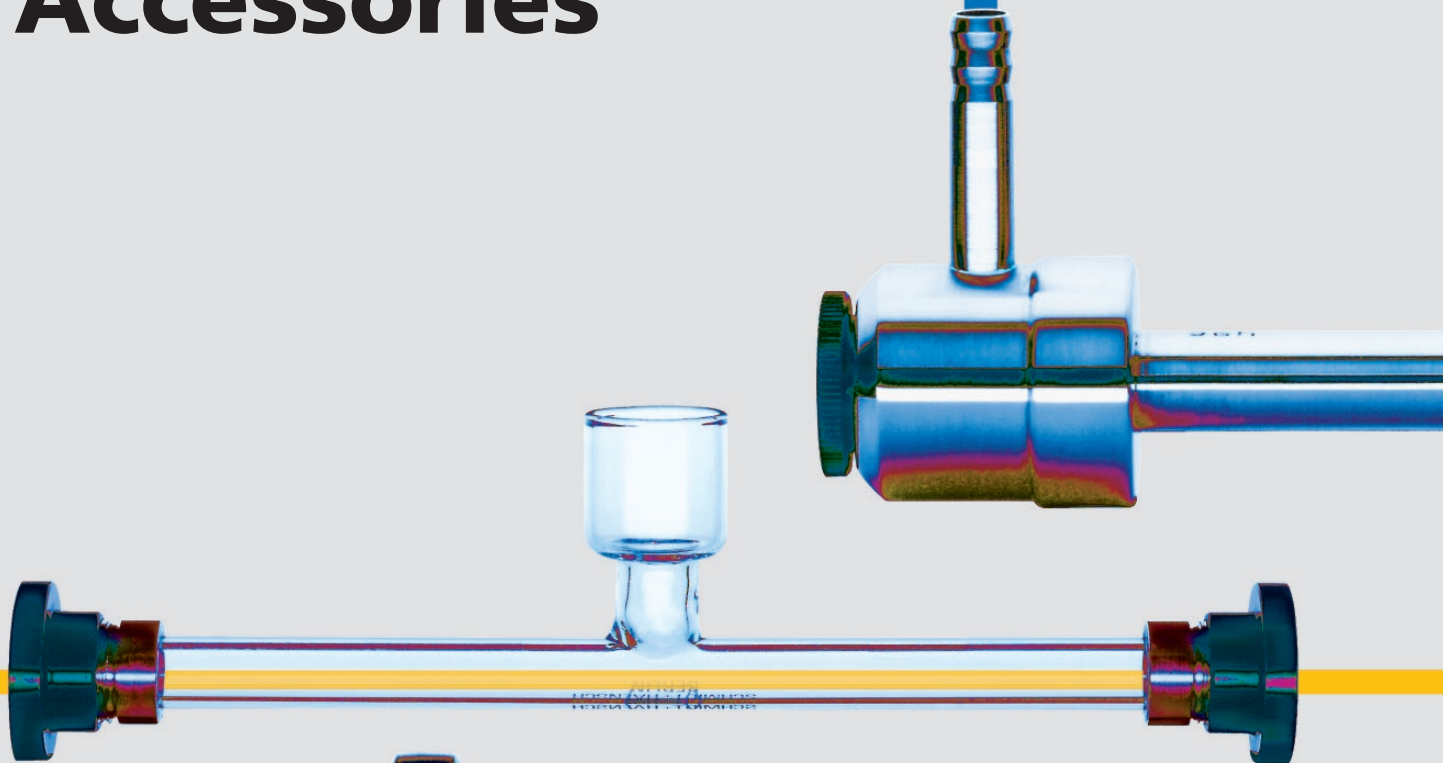

SCHMIDT + HAENSCH

Optisch-elektronische Messinstrumente seit 1864

Polarimeter tubes of SCHMIDT+HAENSCH are designed for different applications. All tubes meet the ICUMSA-recommendations, class A ($\pm 0.01\%$) and the corresponding specifications of the O.I.M.L. and the Australian Standard K 157 regarding of construction and precision.

Hints for the choice of tube length

In general: Longer measuring tubes are better than short ones. Sample heterogeneity and residues from pervious samples have a larger influence on the accuracy of the measurement when applying shorter tubes. Larger sample volumes ensure more stability. Exception: For measuring dark samples it might be necessary to use a shorter measuring tube due to the high absorption of light.

Please consider that SCHMIDT+HAENSCH guarantees for the precision of Polarimeters its only if SCHMIDT+HAENSCH measuring tubes are used.

Flow-through tubes

Flow-through tubes are recommended for testing a large number of samples of the same type.

SCHMIDT+HAENSCH recommends two types of tubes:

- for **manual sample changing** or
- for **automatic sample changing**.

The tubes for manual sample changing work with funnel and riser on the basis of communicating tubes. The sample in the tube is totally displaced by the following sample.

A volume of approx. 60 cm³ is sufficient.

The automatic sample changing in the tubes works using a by-pass or a pump with hose connections. A volume of approx. 60 cm³ is sufficient, too.

Glass sample tubes for single measurements

Glass tube without bubble trap



Length in mm	Volume in cm ³	ID-N°
100	5,5	02595
200	10,5	02596

Glass tubes with bubble trap

Filling through narrow end and short turnover after sealing prevents influence of air bubbles on measurement.



Length in mm	Volume in cm ³	ID-N°
50	8,0	02604
100	10,5	02600
200	16,5	02601

Glass tube with center filling cup

The filling cup allows easy filling and emptying of the tube. Air bubbles escape through the filler opening when moved slightly.




Length in mm	Volume in cm ³	ID-N°
50	3,0	02611
100	6,0	02607
200	12,0	02608

For medical applications:

Length in mm	Volume in cm ³	ID-N°
95,04	5,5	02612
190,09	11,0	02613


New: Temperature sensor Accessories (s.P.4)

Quartz glass flow-through sample tube, thermostatable

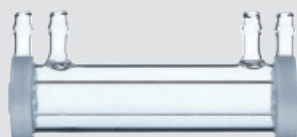

Length in mm	Volume in cm ³	ID-N°
100	7,0	08866

Quartz glass micro tube, thermostatable

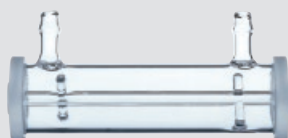

Length in mm	Volume in cm ³	ID-N°
100	1,0	08962

Precision quartz glass flow cell for small sample volume

With adapter for S+H Polarimeter (without image)

Length in mm	Volume in cm ³	ID-N°
10	0,7	09521

► Flow-through sample tubes (stainless steel)

Flow-through tube with funnel and riser,
tube for manual sample changing, not thermostatable



Length in mm	Volume in cm ³	ID-N°
100	12	09212
200	17	09213

Flow-through tube with funnel and riser,
with 6 mm hose connection for water circulation



Length in mm	Volume in cm ³	ID-N°
100	12	09216
200	17	09217

Flow-through tube with funnel and riser,
with integrated temperature sensor, thermostatable



Length in mm	Volume in cm ³	ID-N°
100	12	03820
200	17	03830

Only for the new Polartronic M and H-series and new Saccharomat-models with plug connection inside of the sample compartment. Automatic temperature compensation. **Temperature measuring of the sample**, not of the sample room!

Flow-through tube with hose connection,
automatic sample changing



Length in mm	Volume in cm ³	ID-N°
50	10	09628
100	12	09208
200	17	09209

Flow-through tube, thermostatable version with 6 mm
hose connection for water circulation



Length in mm	Volume in cm ³	ID-N°
100	12	09214
200	17	09215

Micro-flow-through sample tube with 4 mm
hose connection, not thermostatable



Length in mm	Volume in cm ³	ID-N°
10	1,5	02748

► Stainless steel sample tube with center filling cup,
thermostatable



Length in mm	Volume in cm ³	ID-N°
100	12,0	02614
200	17,0	02619

▶ Quartz control plates for calibration of polarimeter (Control of the instrument according to GLP)

Quartz control plates are measuring standards, used as testing and adjusting kits for polarimeters. Thin quartz plates (0.4 to 1.6 mm thickness) are polished with highest precision vertically to the axis and plane parallelly. For reasons of mechanical stability, plates with a low angle of rotation (between -24°Z and $+24^{\circ}\text{Z}$) are produced as double quartz plates (combination of a dextro and a laevorotatory plate).

For each SCHMIDT+HAENSCH quartz control plate a **factory certificate** is supplied, indicating the absolute angle of rotation in angular° and $^{\circ}\text{Z}$ for the wavelength 546, 587, 589 and 882 nm.

The angle of rotation of the quartz plate is directly correlated to its thickness and will be indicated at a temperature of 20°C for the wavelength 589 nm.

Since the angle of rotation is highly influenced by the wavelength while the sugar degrees are not, the following table shows the correlation between the different values of the certificate so that the calculation of the angle of rotation of a quartz plate is possible from the sugar degree indicated below for different wavelengths at 20°C .

Wavelength in nm	546.23	589.44	882.6
Conversion factor from $^{\circ}\text{angular}$ to $^{\circ}\text{sugar}$	2.46	2.89	6.76
Conversion factor from $^{\circ}\text{sugar}$ to $^{\circ}\text{angular}$	0.41	0.35	0.15

The required (or needed) value must be specified in your order!

Single quartz control plate	ID-N°
Standard sugar degree: $-30, +25, +45, +50, +75, +95, +99^{\circ}\text{Z}$ ($\pm 1^{\circ}\text{Z}$)	01817
Sugar degree in range of: -99 to -25 and $+25$ to $+99$ ($\pm 1^{\circ}\text{Z}$)	03227

Double quartz control plate	ID-N°
Standard sugar degree: $+15$ and $+20^{\circ}\text{Z}$ ($\pm 1^{\circ}\text{Z}$)	01818
Sugar degree in range of: -24 to $+24^{\circ}\text{Z}$ ($\pm 1^{\circ}\text{Z}$)	03228

▶ Accessories for quartz control plates	ID-N°
Temperature sensor for automatic temperature correction. Only to be used with new M- and H-series and new Saccharomat® models with plug!	07277



Replacement of the factory certificate (only for quartz plates in good shape and no changes are necessary)
Official certificate for quartz control plates can be supplied by the German Institut of Metrology (PTB) on their actual charges

▶ Accessories for quartz glass tubes and flow-through tubes	ID-N°	
Cover glass for polarimeter tubes, diameter 15.5 mm	01811	
Elastic insert for polarimeter tubes, diameter 15.5 mm	01813	
Cover glass for polarimeter tubes, diameter 23.7 mm	01812	
Elastic insert for polarimeter tubes, diameter 23.7 mm	01814	
Cleaning brush for sample tubes	01838	
Thermometer; 0 to 40°C ; 0.1°C grating; for tubes with center filling cup	01816	
Rubber stopper, with hole to insert thermometer	08341	
New: Temperature sensor for automatic temperature correction for glass tubes with center filling cup (also available for old tubes), can be used with all new M- and H-series and new Saccharomat® models		
	for 100 mm tubes	02624
	for 200 mm tubes	02625

SCHMIDT+HAENSCH GmbH & Co.
Waldstraße 80/81
D-13403 Berlin
Germany
Phone: +49 30 / 41 70 72-0
Fax: +49 30 / 41 70 72-99
e-mail: sales@schmidt-haensch.de
www.schmidt-haensch.com



Subject to modification without notice
Änderungen vorbehalten
05/04